




Oppo 4030 knee immobilizer – Wrap around → Open patella → Rigid

DVA RAP code AR02

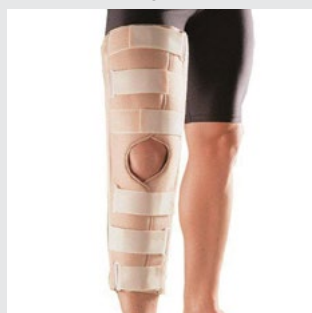

Features:

- provides varying degrees of support and stabilization for the leg and knee.
- foam padding reduces slippage and will not fray.
- aluminum posterior stays retain proper position for comfort.
- hook and loop closure provides quick, secure fit.
- available in three heights, 45.7cm, 50.8cm and 58.4cm

Application:

- post surgical use
- muscular contracture

To Apply:

- wrap immobilizer around the knee.
- position the aluminum stays on each side of knee.
- be sure the opening of support should be on kneecap and secure the hook and loop closures.
- thread through the buckles and secure the upper and lower hook and look straps to desired compression.

46cm length	51cm length	58cm length	size	cm
OPP4030S-46	OPP4030S-51	OPP4030S-58	S	33-38.1
OPP4030M-46	OPP4030M-51	OPP4030M-58	M	38.1-45.7
OPP4030L-46	OPP4030L-51	OPP4030L-58	L	45.7-53.3
OPP4030XL-46	OPP4030XL-51	OPP4030XL-58	XL	53.3-61

Measure the circumference of the thigh approximately 10cm above the knee.



Oppo tri-panel knee immobiliser – Wrap around → Open patella → Rigid

DVA RAP code AR02



Features:

- universal knee immobiliser
- 3 panel design
- 2 side panels velcro on to the middle panel to allow a range of circumferential sizing.
- one size fits all (length: 51cm)

Application:

- post surgical use
- muscular contracture

AUPHYKI	Universal size
---------	----------------

Post-op knee immobilizer – Wrap around → Open knee → Rigid

DVA RAP code AR02



Features:

- removable medial, lateral and posterior aluminium stays can be contoured to the shape of the leg for a comfortable and secure fit
- the 3 panel design (including the popliteal variation) allows adjustment of the brace to cater for all sizes of leg throughout the course of the patient's treatment.
- will fit calf to 46.5cm, and thigh to 56.5cm

Application:

- allows full immobilisation of the knee in full extension, while providing mediolateral support

Product code	size	Length cm
SM800	universal	Short 40
SM813	universal	Med 48
SM814	universal	Long 58