




Oppo 1029 patella tendon strap – Wrap around (black)

DVA RAP code AR02

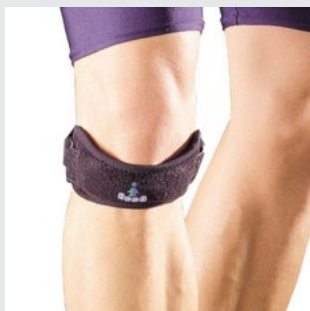

Features:

- extra pad provides compression on target area
- non-slip backing keeps wrap from sliding out of position
- great for use after an injury or to help prevent injuries
- reduces vibrations in the femoral groove

Application:

- patellar tendinitis
- osgood-Schlatter syndrome

To Apply:

- position the pad below the patella
- fasten hook and loop closure strap securely

OPP1029	One size	31.8-40.6
---------	----------	-----------

Measure around the knee joint.



Oppo 1429 jumpers strap w/pad – Wrap around (beige)

DVA RAP code AR02



Features:

- specially shaped silicone insert offers extra pressure on the patella
- adjustable hook and loop strap allows you to apply greater tension if needed
- minimizes potential for future injuries

Application:

- patellar tendinitis
- jumper's knee
- osgood-Schlatter syndrome

To Apply:

- wrap around the knee and position the silicone pad below the knee cap
- fasten with hook and loop closure
- thread through the buckle and secure the elastic hook and loop strap to desired compression

OPP1429	One size	31.8-40.6
---------	----------	-----------

Measure around the knee joint.



McDavid A414 patella knee strap – Wrap around Neoprene (black, beige)

DVA RAP code AR02



Features:

- a semi-tubular buttress in a thermal strap with velcro closure
- for compression and pressure on the patellar tendon
- exerts pressure over the patella tendon
- one size fits most

Application:

- suitable for chondromalacia patella, patellar tendonitis, Osgoode-Schlatter, general aches & pains behind the knee, runner's knee, bowler's knee, & patella bursitis

Black	Beige	Size
NHA414B	NHA414S	Universal size


Oppo 1028 patella tracker support – Wrap around → Open patella

Neopren (black)

DVA RAP code AR02

**Features:**

- prevents patella pain and weakness due to instability
- helps maintain proper patella position
- helps prevent Chondromalacia Patella
- helps prevent patella subluxation and dislocation
- pressure pads and straps help secure the patella in the proper position
- does not interfere with normal motion, even when straps are crossed
- the trapezium helps provide comfort
- the pad design ensures comfortable fit and provides even pressure to the side of the knee

Application:

- patella Instability, tilting, subluxation
- patellofemoral pain symptom
- chondromalacia patella
- post operative treatment

To Apply:

- place the opening over the kneecap so the pad is placed on the outside of the kneecap (In some cases your doctor or the therapist may recommend placing the pad on the inside. Please refer to your doctor or therapist's instructions.)
- attach the short hook and loop straps in the knee fold. Smooth down the material
- put the lower long strap in the upper plastic ring, then pull back and attach the hook and loop straps. (The long straps regulate the pressure from the pad against the kneecap.)
- the pad should now feel tight against the kneecap

Product Code	Size	cm
OPP1028XS	XS	30.5 - 33.0
OPP1028S	S	33.0 - 35.6
OPP1028M	M	35.6 - 38.1
OPP1028L	L	38.1 - 40.6
OPP1028XL	XL	40.6 - 43.2
OPP1028XXL	XXL	43.2 - 45.5

Measure around the mid-section of the knee, while in a standing position.


Body Assist PSB patella stabilizing brace – Wrap around → Open patella

Neoprene (black)

DVA RAP code AR02

**Features:**

- designed according to knee cap taping technique, locking the knee cap
- anatomically shaped to provide orthopaedic support
- retains some body warmth & can increase blood circulation
- a great brace for anyone who needs to tape their knee cap

Application:

- for general deficiency in knee cap tracking movement & assists in the treatment of chondromalacia patella, general aches & pains behind the knee, runner's knee, bowler's knee, & patella bursitis

Product Code	Size	cm
NHPSBCHILD	Child	to 32
NHPSB	Regular	32-42
NHPSBXLG	XLarge	42 +

Measure around the knee joint.


Rehband 8385 genu therma tendon – Wrap around

Neoprene (grey)

DVA RAP code AR02

**Features:**

- relieves the pressure on the patellar tendon insertion
- provides pain relief
- knee strap with a soft pad inside
- excellent for sports and exercises when you have pain behind or below the kneecap
- can be used for left or right

Application:

- mb. Schlatter
- jumper's knee
- anterior knee pain

Product Code	Size	15 cm below mid patella	15cm above mid-patella
REH8385XXS	XXS	29-32	36-40
REH8385XS	XS	32-35	40-44
REH8385S	S	35-38	44-48
REH8385M	M	38-41	48-52
REH8385L	L	41-44	52-56
REH8385XL	XL	44-48	56-61
REH8385XXL	XXL	48-51	61-67

Measure around leg, 15cm above, and 15cm below middle of knee.



Rehband 8386 genu therma chondro – Wrap around → Open patella

Neoprene (grey)



Features:

- relieves the pressure on the patellar tendon insertion
- patella tendon strap with an anchor
- the adjustable anchor above the knee prevents the brace from migrating or rotating
- excellent for sports and exercises when you have pain behind or below the kneecap
- can be used for left or right

Application:

- mb. Schlatter
- jumper's knee
- anterior knee pain

DVA RAP code

AR02

Product Code	Size	15 cm below mid patella	15cm above mid-patella
REH8386XXS	XXS	29-32	36-40
REH8386XS	XS	32-35	40-44
REH8386S	S	35-38	44-48
REH8386M	M	38-41	48-52
REH8386L	L	41-44	52-56
REH8386XL	XL	44-48	56-61
REH8386XXL	XXL	48-51	61-67

Measure around leg, 15cm above, and 15cm below middle of knee.



Rehband 8360N genu carezza patella stabilizer – Wrap around → Open patella

Elastic (grey)



Features:

- keeps patella in correct position
- prevents patella from luxation
- supports sensory motor control
- flexible patella stabilizer
- c-shaped pad, which is adjusted by two separate straps
- innovative material that combines compression and breathability
- comfortable to wear
- excellent fit underneath clothes
- can be used for left or right

Application:

- patellofemoral pain syndrome (e. g. chondromalacia, chondropathy of the patella, beginning osteoarthritis of the patellofemoral joint)

- patellofemoral malalignment
- stabilization after patella luxation/subluxation
- stabilization after lateral release surgery

DVA RAP code

AR02

Product Code	Size	15 cm below mid patella	15cm above mid-patella
REH8360NXXS	XXS	29-32	36-40
REH8360NXS	XS	32-35	40-44
REH8360NS	S	35-38	44-48
REH8360NM	M	38-41	48-52
REH8360NL	L	41-44	52-56
REH8360NXL	XL	44-48	56-61
REH8360NXXL	XXL	48-51	61-67

Measure around leg, 15cm above, and 15cm below middle of knee.



Oppo 1033 knee support patella stabilizer – Pull on → Open patella → Stabilising straps

Neoprene (beige)



Features:

- neoprene circle pad around the patella provides protection and stabilization to prevent patella displacement.
- four neoprene straps provide support directly to patella area.
- two medial and lateral spiral stays provide additional support and stability for the knee.
- helps relieve pain from tendinitis and arthritis.
- retains body heat to warm knee joint and increases blood circulation.

Application:

- moderate arthritis
- knee strain/ sprain
- knee instability
- patellar instability

To Apply:

- hold the top of support using both hands, one on each side.
- insert foot into knee stabilizer and pull up to the knee.
- position the opening over the kneecap.
- be sure the two medial and lateral stays are on each side of knee.
- pull the medial straps above and below the patella of one side to tighten and secure them to the lateral side of the knee to desired compression.
- repeat step 5 to secure another side of straps.

DVA RAP code

AR02

Product Code	Size	cm
OPP1033S	S	33.7-37.5
OPP1033M	M	37.5-40.6
OPP1033L	L	40.6-45.1
OPP1033XL	XL	45.1-50.2

Measure around the mid-section of the knee, while in a standing position


Rehband 8387 genu therma combi – Pull on → Open Patella

Neoprene (grey)

DVA RAP code AR02

**Features:**

- relieves the pressure on the patellar tendon insertion
- adjustable patella tendon strap
- the adjustable anchor above the knee prevents the brace from migrating or rotating
- excellent for sports and exercises when you have pain behind or below the kneecap
- can be used for left or right

Application:

- mb. Schlatter
- jumper's knee
- anterior knee pain

Left	Size	15 cm below mid patella	15 cm above mid patella
REH8387XXS	XXS	29-32	36-40
REH8387XS	XS	32-35	40-44
REH8387S	S	35-38	44-48
REH8387M	M	38-41	48-52
REH8387L	L	41-44	52-56
REH8387XL	XL	44-48	56-61
REH8387XXL	XXL	48-51	61-67

Measure around leg, 15cm above, and 15cm below middle of knee.