



SWELLSPOTS

MANUFACTURED BY SOLARIS

Scar and Edema Pads that
provide simple solutions.

Swell Spots

Quick Reference Guide

Solaris Quick Reference Guides are your helpful resource for basic questions and information about our various product lines. If you're not able to find what you're looking for or would like more information, please call us at 855.892.4140 or e-mail us at info@solarismed.com.

TABLE OF CONTENTS

- 2 Contact Information
- 3 About Swell Spots
- 4 Swell Spot Details & Diagrams

CUSTOMER SERVICE

Customer Service Representatives are available to assist you from 7 AM to 7 PM CST Monday through Friday.

855.892.4140 (Toll-free)

414.918.9189 (Fax)

info@solarismed.com

Solaris, Inc.

6737 West Washington St., Suite 3260

West Allis, WI 53214 USA



HIGH PERFORMANCE FABRIC

Our grey and black fabrics promote comfort with their breathable, wickable, and thermal release properties. The grey side has tiny microfibers that grip skin and encourage gentle tissue stretch to help open the initial lymphatics, ensuring comfort with therapy.



ENDLESS OPTIONS

Swell Spots can be placed almost anywhere on the body and tucked into bandage liners and day or night-wear garments for specific attention to problematic areas. There's no limit as to how or where they can be used.



FOAM CHIP TECHNOLOGY

Unique foam particles help create localized tissue pressure differentials that stimulates interstitial fluid movement, prevents stagnation, and softens indurated tissue.



CARE INSTRUCTIONS

Machine wash and dry using mild detergent.

MATERIALS

Latex-free fabric and foam.

INDICATIONS

- Lymphedema
- Other forms of edema


CONTRAINDICATIONS

- Arterial Insufficiency or degeneration
- Deep Vein Thrombosis
- Untreated or congestive heart failure
- Localized or systemic infection
- Untreated Cancer
- Absent or severely impaired sensation

CURVED (SP-01)

Measures approximately
5 cm x 12 cm


Primarily used for hard to
accommodate areas (ankles
and knee), but also effective
around the elbow, inferior breast
scarring, around the eyes, or
behind the ear.

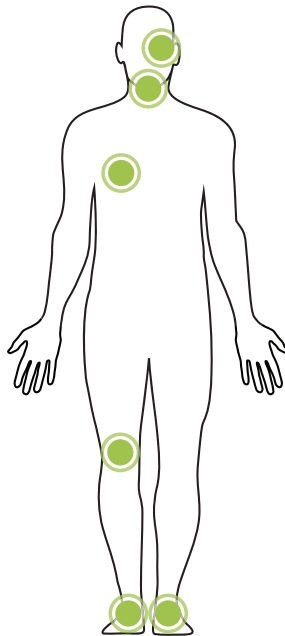

DORSAL (SP-02)

Available in Small and Large.
Small (SP-02-S) : 9 cm x 14.5 cm
Large (SP-02-L) : 9 cm x 18 cm



Designed for use on the
back of the hand or foot.

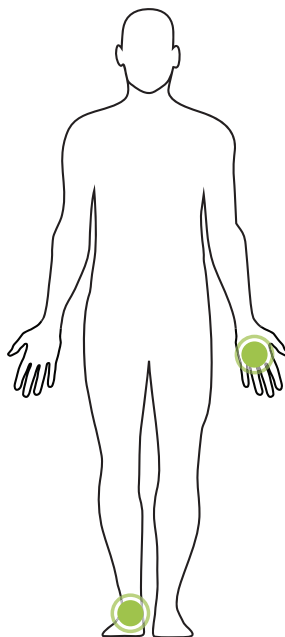

PALMAR (SP-03)

Available in Small and Large
Small (SP-03-S) : 7 cm x 8.5 cm
Large (SP-03-L) : 8 cm x 9.5 cm


Primarily for use the palm; also effective on the back of the hand, and lumpectomy scarring when used in a bra.

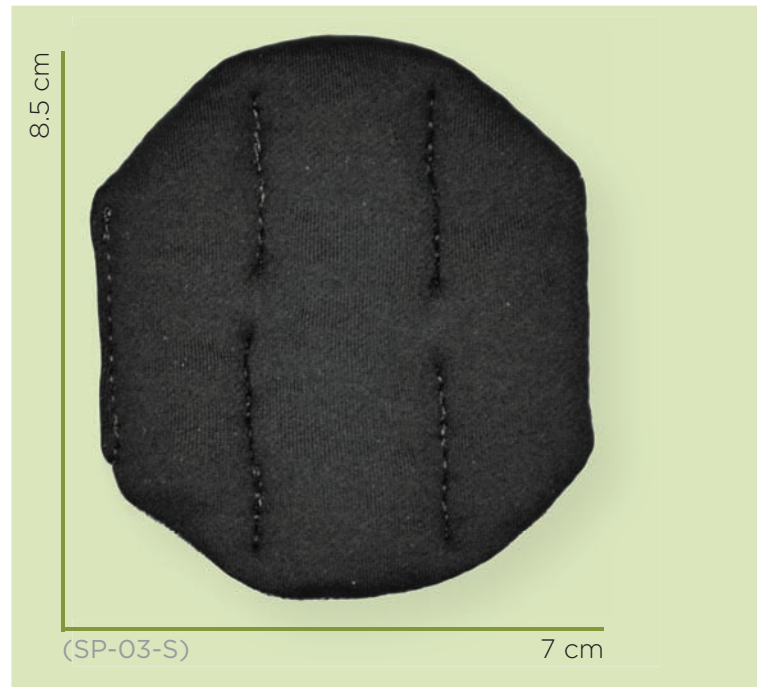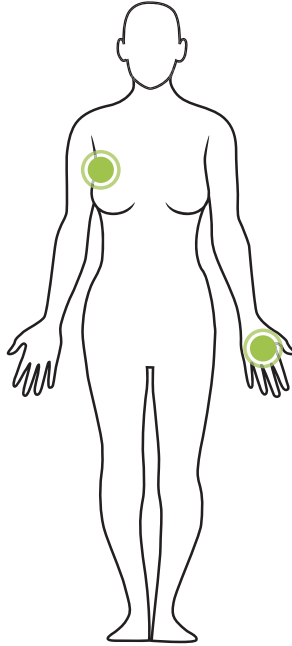

OVAL (SP-04)

Measures approximately
8 cm x 19 cm


Versatile and easy to use with a simple, multi-purpose design allowing for placement in many areas, such as inside hand splints, along the ulnar aspect of the forearm, and labia swelling.

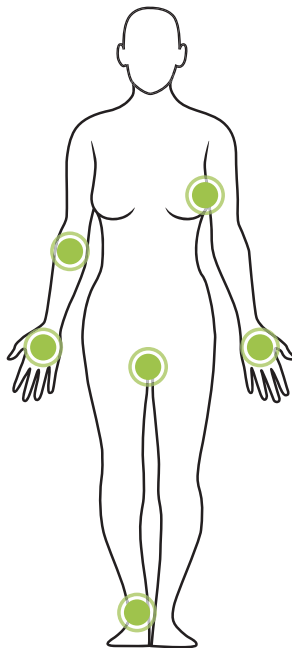

KATINA (SP-05)

Available in Small and Large.
Small (SP-05-S) : 10 cm x 38 cm
Large (SP-05-L) : 18 cm x 42 cm

Tapered shape useful in treating difficult areas, such as ulnar aspect of forearm, IT band of lateral thigh, and anterior surface of calf.

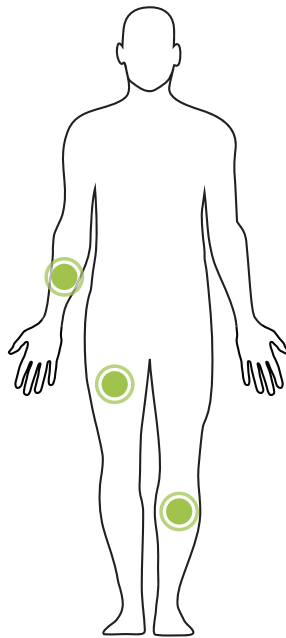


SOBOL (SP-06)

Small (SP-06-S) : 18 cm diameter
Large (SP-06-L) : 23 cm diameter

Addresses post-mastectomy fibrotic tissue on chest wall due to radiation and/or surgical incision scarring.


Also useful for hip, anterior or posterior knee, and abdomen when used under bandages and compression garments.

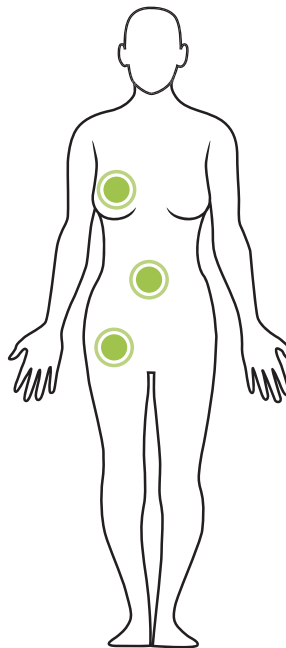

JOJU (SP-07)

Measures approximately
5 cm x 8 cm



Specifically addresses swelling
of thenar web space; also
effective at softening scars
from carpal tunnel surgery.

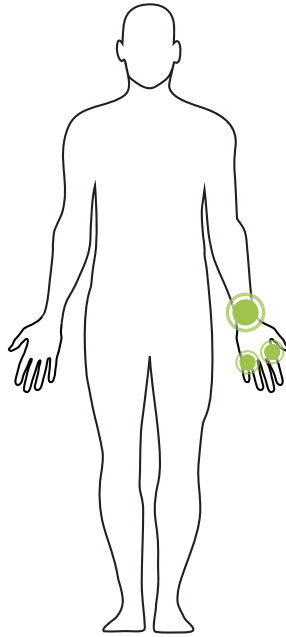

KIMBE (SP-08)

Measures approximately
15 cm x 25 cm

Appropriate for reducing
swelling over the shoulder
blade, as well as over the
groin and neck area.

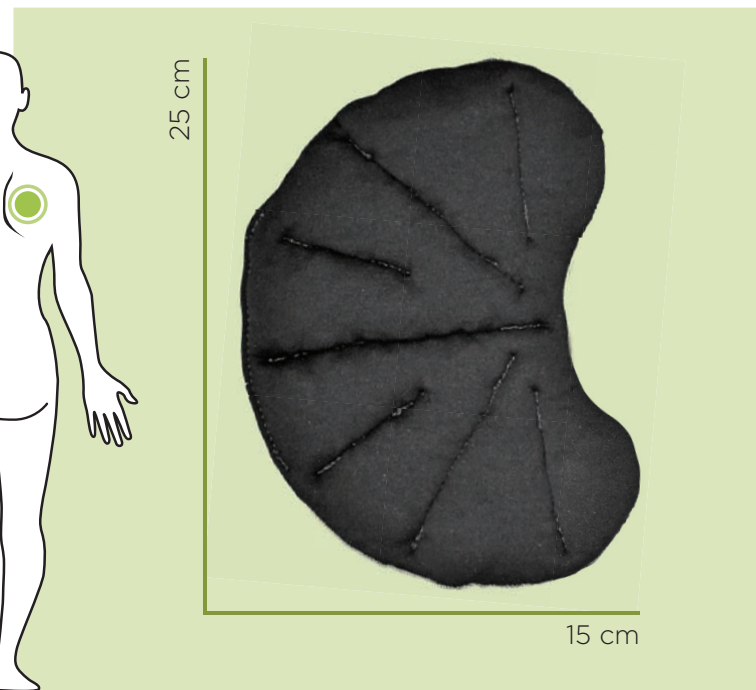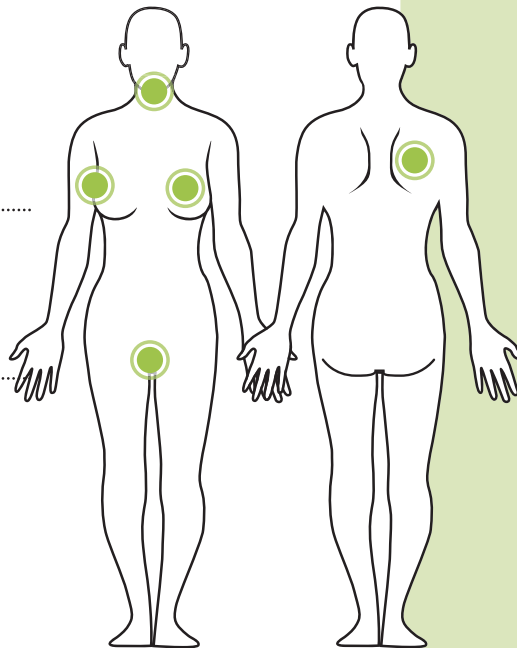


BREAST (SP-09)

Small (SP-09-S) : A-B cup
Medium (SP-09-M) : C-D cup
Large (SP-09-L) : D-DD cup
*Sizes are approximate

Designed specifically for swelling and fibrosis of breast area after lumpectomy or reconstructive surgery.


Seamless design allows spot to be used on either left or right sides.

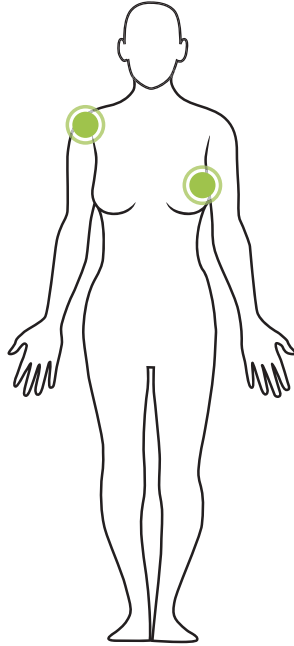

VERSI-FACE (SP-10)

Measures approximately
17 cm x 26 cm

Designed to address swelling around eyes; also effective when placed under chin covering the neck, or rotated to cover mandible and cheeks.

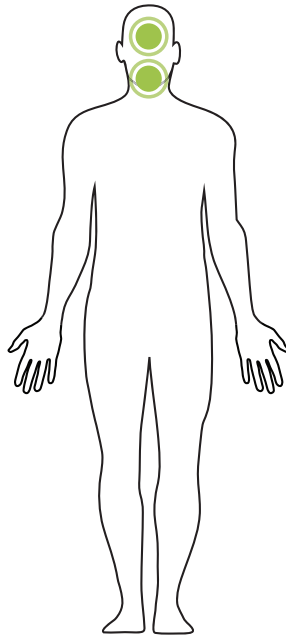


THORA-THIGH (SP-11)

Measures approximately
22 cm x 50 cm

Designed for asymmetrical chest
coverage, groin, medial thigh, and
ankle areas.

Effective on radiation fibrosis
and dispersing ending pressure
of bandages and leg garments.

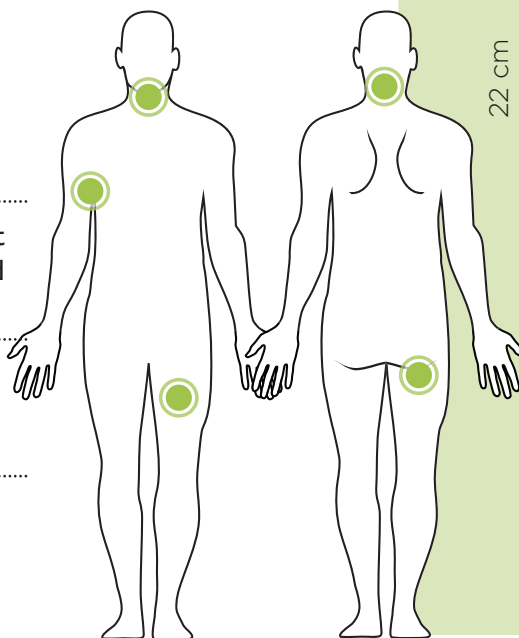

22 cm




50 cm

REYNOLDS (SP-12)

Pad measures approximately
7.5 cm x 16.5 cm

Reynold's strap constructed
with Velcro-receptive material
that hooks to Spot and can be
placed in multiple locations.

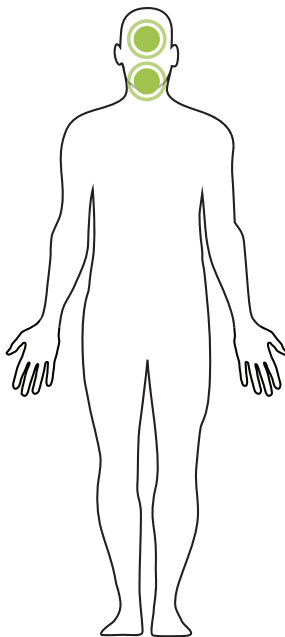

CRESCENT (SP-13)



Available in Small, Medium, and Large.

Small (SP-13-S) : 54 cm L x 6.5 cm H x 9 cm D

Medium (SP-13-M) : 60 cm L x 6.5 cm H x 11 cm D

Large (SP-13-L) : 71 cm L x 6.5 cm H x 12 cm D

Effective when placed in skin folds at reducing maceration and odor, softening upper/lower lobules, and supporting lax abdominal tissue.

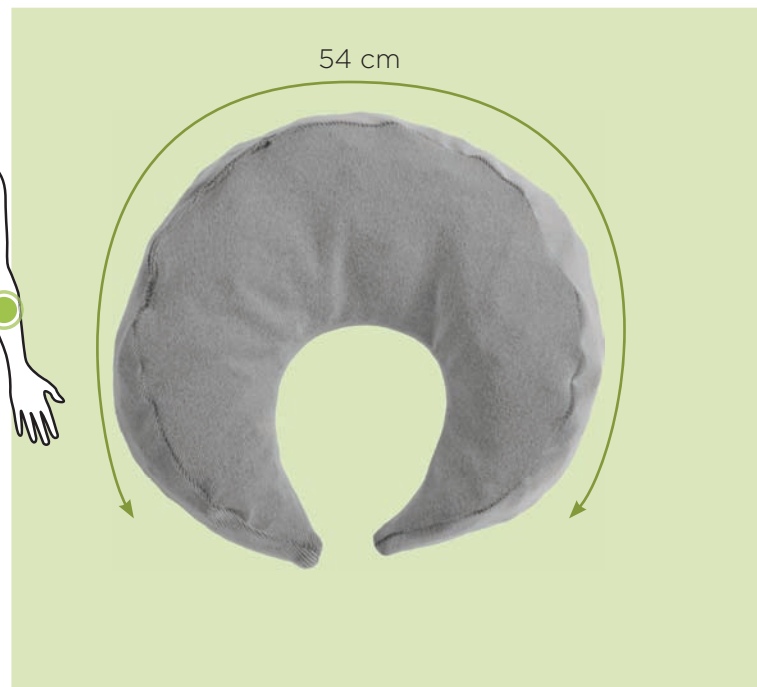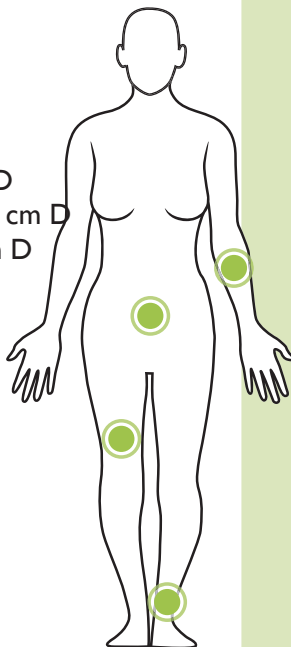

WE'A TOE (SP-14)

Available in Small, Medium, and Large.


Small (SP-14-S) : 22 cm MTP circumference

Medium (SP-14-M) : 25 cm MTP circumference

Large (SP-14-L) : 32 cm MTP circumference

Exclusive digit spacers offer comfort and compression for toes individually.


Must be able to reach toes to don.

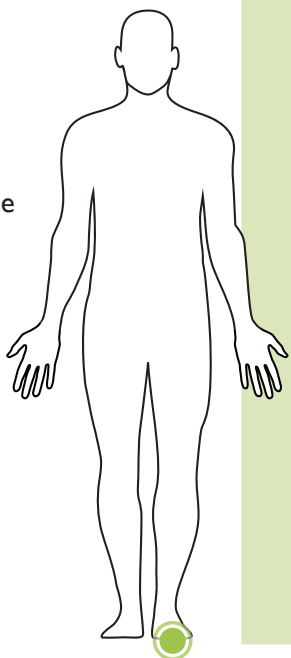

LABIA (SP-15-F)

Measures approximately
30 cm x 50 cm


Labia Swell Spots feature tapered design to enhance comfort while applying gentle compression and lower abdominal coverage to direct fluid away from genital area; can be worn inside or over underwear.

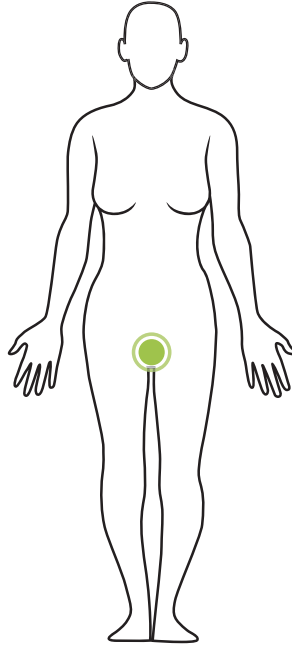

SCROTUM (SP-15-M)

Measures approximately
30 cm x 50 cm


Scrotum Swell Spot designed to gently cradle scrotum; accommodates up to 27 cm of scrotal circumference, though further customizing is available.

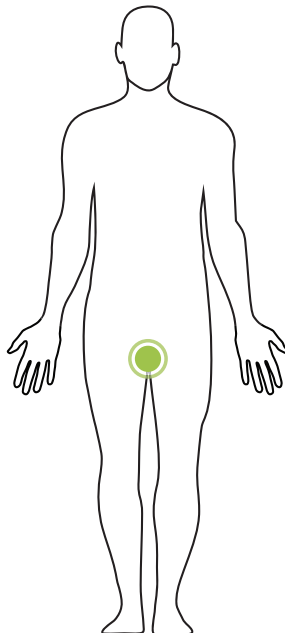

LATERAL BRA (SP-16-A)

Lateral Bra (SP-16-A) measures approximately 10 cm x 16 cm


Designed to tuck inside compression bra pockets or regular bra to reduce scarring and protect delicate areas from bra chest strap pressures.

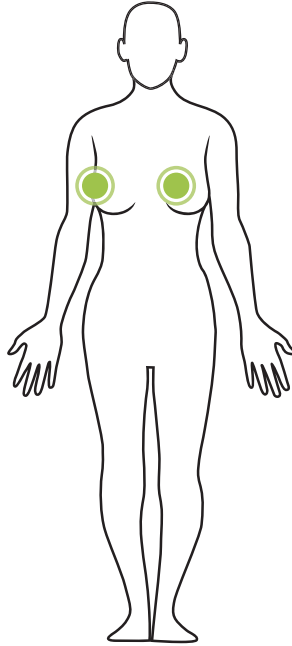

FULL BRA (SP-16-B)

Available in Small, Medium, and Large.
Small (SP-16-B-S) : A-B Cup
Medium (SP-16-B-M) : B-C Cup
Large (SP-16-B-L) : C-D Cup

Designed for the Belisse bra or to tuck into sports bras with front openings.


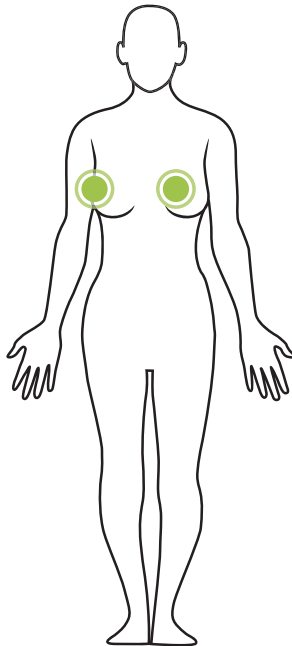


LT KNEE SHAPER (SP-17)

Patella girth range 35 cm to 65 cm

Helps contour/shape knee while protecting delicate structures.

Use for knees or large elbows to pad bony prominences, diffuse bandaging pressures, and tissue induration.

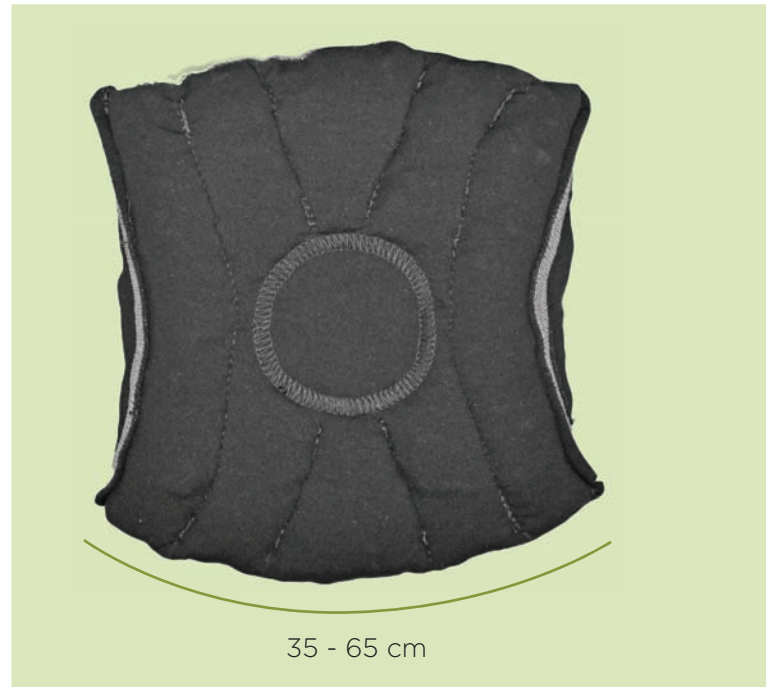
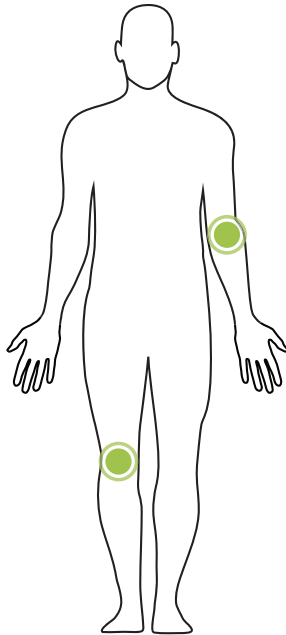


RING (SP-18)

Measures approximately
54 cm L x 6.5 cm H x 9 cm D

Available with or without
hook and loop closure.

Designed to fill crevices at the
ankle and wrist; also effective
at treating post-surgical
scarring from hysterectomies.

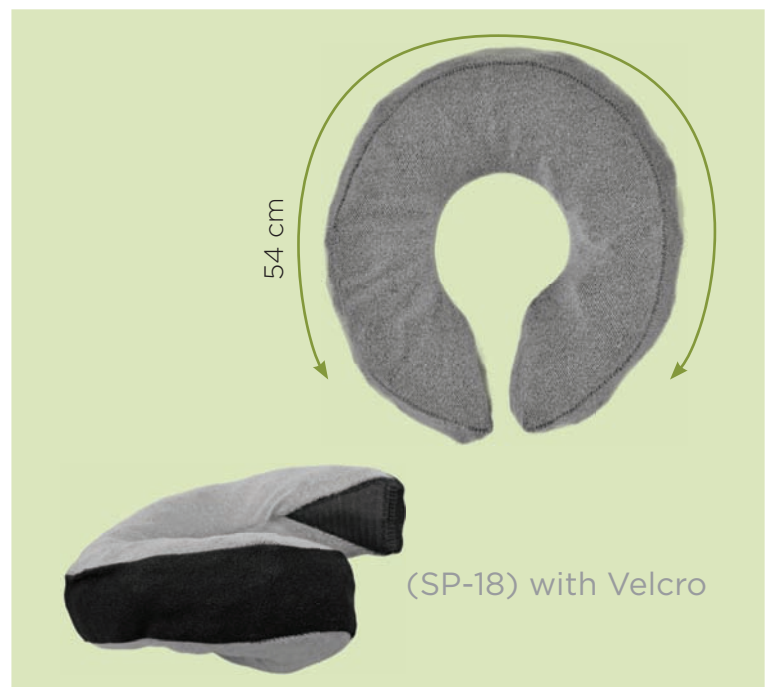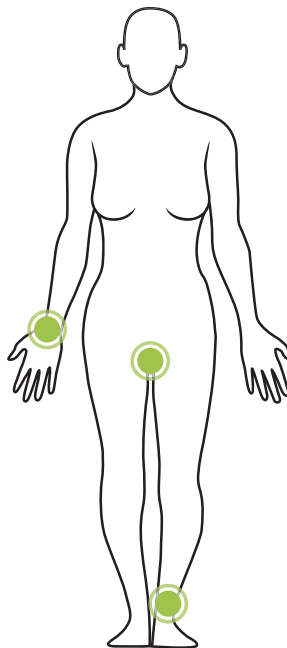


SWEET SPOT (SP-19)

Measures approximately
7 cm x 14 cm

Features one-sided non-skid
surface to prevent repositioning
during normal movement.

Mimics size and shape of normal
panty liner to be worn inside,
over, or between underwear
for maximum comfort.


Addresses residual and transient
pelvic floor swelling in pregnant
and postpartum women, as well
as perineum and labia swelling.

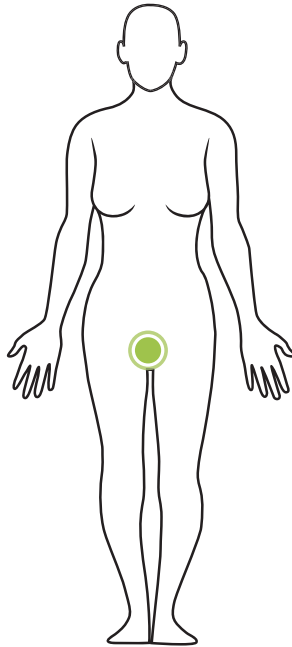

TIC-TAC (SP-20)

Small (SP-20-S) : 18 cm x 18 cm
Medium (SP-20-M) : 28 cm x 28 cm
Large (SP-20-L) : 38 cm x 38 cm

Works to address larger areas of fibrosis.

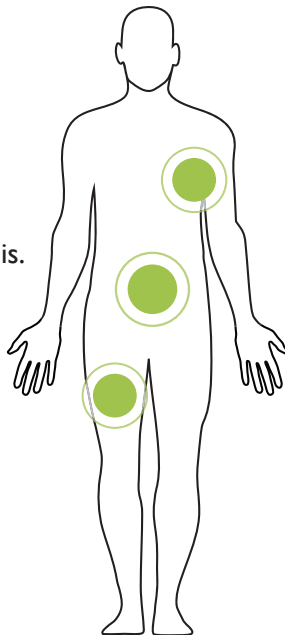



Connect with us
by liking our Facebook page!



facebook.com/solarisinc



twitter.com/solarisinc



youtube.com/solarisinc

Solaris, Inc.
6737 W. Washington Street, Suite 3260
West Allis, WI 53214
855.892.4140 (Toll-free)
414.918.9189 (Fax)
info@solarismed.com